

Plinth 2000

Model 93CDT Tilting Divided Leg Podiatry Chair

Overview:

Two motor podiatry treatment chair with a full 90° leg drop, shorter seat and trendelenburg tilt greatly enhancing patient comfort and access whilst doing heel work. This model features smooth gas strut operated back rest and divided & extending leg rests. The stay level, swing back arm rests are removable to allow multi-disciplinary use of the chair in shared clinics. A wide range of options can be added to expand the use of this chair.


Standard Features

- UK Manufactured
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Electric operation with trendelenburg tilt
- Comfort foam
- Easy operation
- Low running costs
- Hardwearing anti-microbial vinyl
- Safe working load 250kg
- Retractable castors
- Removable fold back arm rests
- Divided & extending leg rests
- Hoist clearance 130mm

Overall Dimensions

Length: 176cm (69") approx
 Width: 63cm (25") approx
 Height: 52cm - 101cm approx (20" - 40") approx

Upholstery Dimensions

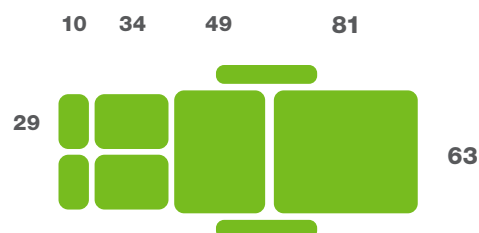

Recommended Options:



Double Foam



Foot Switch



Battery Back-up



D Cushion



Motor Back




Paper Roll Holder

NHS Approved Supplier




t: +44 (0)1449 767887
 e: sales@plinth2000.com
 www.plinth2000.com

Operational Parameters

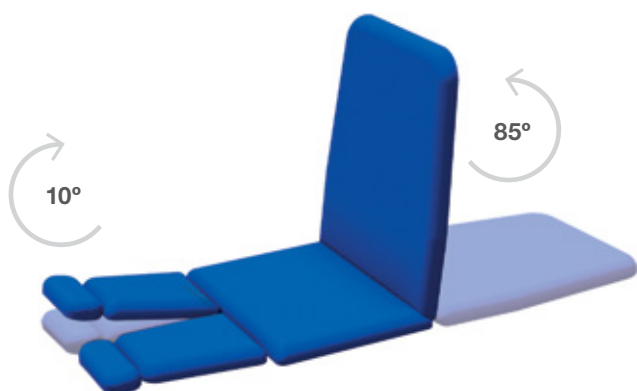

Trendelenburg Tilt 20°

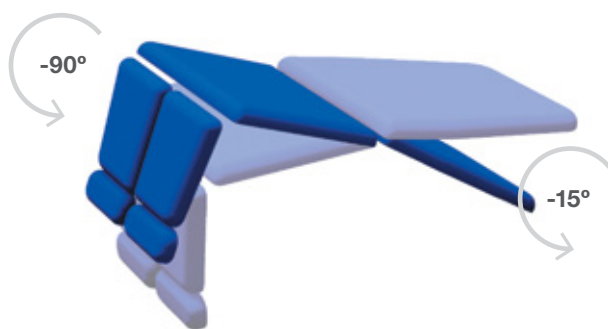

Model 93CDT - Electric version

For Your Peace of Mind

- Lifetime frame & 5 year parts guarantee (accidental damage to upholstery excluded)
- Optional extended warranty to 10 years
- Average lifespan of 16 years
- Designed and manufactured in-house in the United Kingdom
- CE marked
- Externally audited to ISO9001:2008 by LRQA

Standards

- IP Rating - IPx4
- EEC Classification - MDD93/42/EEC Annex IX Class 1
- Electrical Standard - IEC601 Part 1 + 60601. P (EN)
- Vinyl Flame Retardant To CRIB 5 (BS7176 Including Foam)
- Foam To BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]

Delivery Promise

- Delivery to UK mainland in just 7-10 working days
- Faster service maybe available if required

Upholstery Vinyl

- Flame retardant to CRIB 5
- Patented anti-microbial technology proven to inhibit the growth of E-Coli and MRSA
- Excellent levels of wear and longevity
- Fully washable and impermeable, to aid infection control
- Suitable for cleaning and disinfecting with dilute solutions of Hypochlorite (10,000ppm)
- Please note: Unless specified, beige will be supplied.



Plinth 2000 Ltd.
Wetheringsett Manor
Wetheringsett
Stowmarket
Suffolk, IP14 5PP

t: +44 (0)1449 767887
f: +44 (0)1449 766122
e: sales@plinth2000.com
www.plinth2000.com

Plinth | 2000