



Overview:

New model with twin column lift, allowing greater strength and optimum positioning of bariatric patients for a variety of gynaecological and colposcopy procedures. The reverse Trendelenburg feature greatly aids patients in accessing the chair whilst full electrical operation of height, backrest and Trendelenburg tilt increase patient comfort, positioning and clinical access. Supplied as standard with powered stirrups, removable arm supports, head cushion and stainless steel tray.

Benefits


- Lifetime frame & 5 year parts guarantee - see website for further details
- Electric operation of height, backrest and 2 way tilt to +/- 20°
- Seat section with cut-out area and option of extension for patient positioning and recovery
- Motorised backrest for supporting the patient in various seated positions
- Safe working load of 320kgs (704lb / 50 stone)
- Large 125mm twin castors with central locking and steering facility
- Supplied with powered stirrups to assist patient positioning
- Battery backup included incase of mains power failure
- Low patient transfer height of 54cm (21")
- Wide choice of options and accessories to meet user requirements.




Overall Dimensions

Length: 143cm (56") approx
 Width: 78cm (31") approx
 Height: 54cm - 94cm approx
 (21" - 37") approx

Upholstery Dimensions

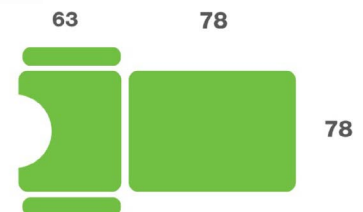


Options & Accessories:

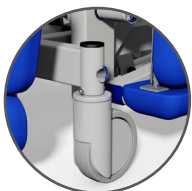

Electronic Braking Castors



Foot Switch

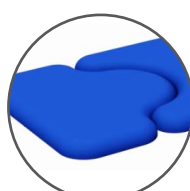

Seat Extension



Matching Operators Chair

NHS Approved Supplier




t: +44 (0)1449 767887
 e: sales@plinthmedical.com
 www.plinthmedical.com

Operational Parameters



Trendelenburg Tilt 20°

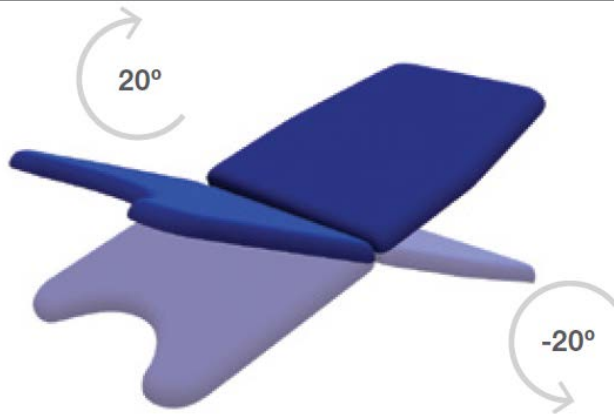

Model 51G - Electric Version

Standard Features

- Designed and manufactured in the UK
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Trendelenburg & reverse Trendelenburg tilt
- Designed and manufactured in the UK
- Vertical dual column lifting operation
- Motorised tilt, height and backrest
- Safe working load and lifting capacity of 320kg
- Large 125mm twin castors with central locking and steering
- Hand switch operated electronic CPR
- Powered stirrups
- Battery back-up
- Low running costs
- Hardwearing anti-microbial vinyl
- Hoist clearance 130mm
- Stainless steel tray, head cushion and folding arm supports included as standard.

Technical Standards

- IP Rating - IPx4
- EEC Classification - MDD93/42/EEC Annex IX Class 1
- Electrical Standard - IEC601 Part 1 + 60601 P (EN)
- Vinyl Flame Retardant to CRIB 5 (BS 7176)
- Foam to BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]
- CE Marked
- Externally audited to ISO9001:2015 by LRQA

Delivery Promise

- Delivery to UK mainland in just 7-10 working days
- Faster service may be available if required
- Overseas delivery 7-10 days plus shipping time

Upholstery Vinyl

- Flame retardant to BS 5852 Ignition Source 5 - CRIB 5 (BS 7176) Medium Hazard
- Independently tested to verify effectiveness at reducing the spread of MRSA, E.coli & other similar bacteria & fungi
- Rated at 500,000 Martindale rubs without damaging or marking the top layer, making it suitable for all high-traffic, heavy-duty areas
- Phthalate free
- Please note: If the upholstery colour is not specified, 'Almond' will be supplied



Jasmine



Battleship



Cobalt



Cocoa



Grape



Mulled Wine



Pillarbox



Gingersnap



Candy



Applemint



Sapphire



Atlantic Blue



Lupin



Denim



Cool Blue



Rainforest



Wasabi



Citrus Green



Marigold



Almond